



## PowerGrip Bipolar Laparoscopic Instruments

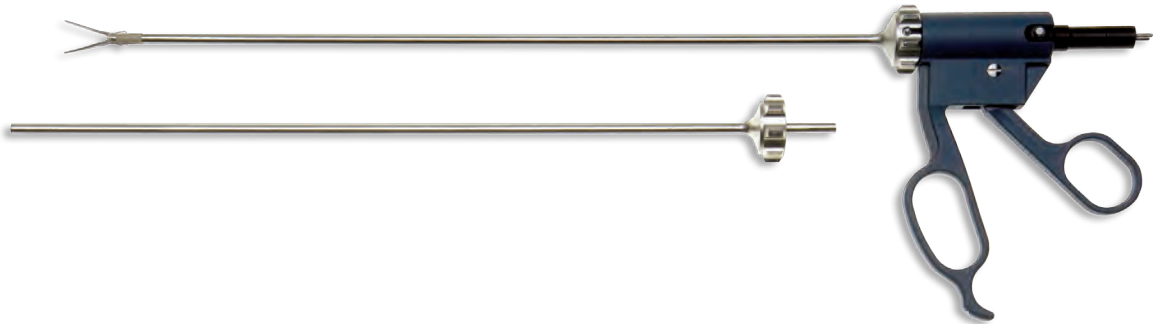


### PRODUCT HIGHLIGHTS

- PowerGrip bipolar laparoscopic instruments facilitate precise dissecting, grasping, cutting and coagulation.
- The PowerGrip bipolar handles and inserts are reusable.
- Jaw mechanism enables adjustable opening and closing of the jaws with very high pressure while grasping and cutting.
- Smooth handle operation provides for an optimal tactile feel and mechanical advantage.
- Instrument jaws are insulated to the end of the grasping and cutting zone in order to avoid unintentional coagulation.
- Rotation knob allows electrode insert to be rotated up to 360° with the forefinger.
- Electrode inserts can be interchanged in the sterile field.
- Modular instruments facilitate ease of assembly and cleaning.
- Bipolar instrument cord connection is integral to handle and standardized to prevent improper connections and to eliminate the purchase of additional custom cables.
- Bipolar instrument cord options allow for fast, easy connection to most commonly used ESUs.


## Reusable PowerGrip Handle Assembly

Handle Assembly	Catalog #	
	34cm	45cm
 <p>PowerGrip Handle, 2-Pin</p> <p>BP7500</p>		
 <p>PowerGrip Shaft</p>	BP7550	BP7550-45

## Reusable PowerGrip Bipolar Inserts

Tip Style	Catalog #	
	34cm	45cm
 <p>Bipolar Forceps</p>	BP7601	BP7601-45

## PowerGrip Bipolar Instrument Cords

Cords	Catalog #
 <p>Bipolar Instrument Cord, 2-Pin Angled, Disposable</p>	BP4510 (50 per case)